

## PRODUCT DATA SHEET

Bioworld Technology,Inc.

# NF-kappaB p65 monoclonal antibody

Catalog: MB65963 Host: Mouse Reactivity: Human, Mouse, Rat

#### **BackGround:**

Transcription factors of the nuclear factor (NF-κB)/Rel family play a pivotal role in inflammatory and immune responses. There are five family members in mammals: RelA, c-Rel, RelB, NF-κB1 (p105/p50), and NF-kB2 (p100/p52). Both p105 and p100 are proteolytically processed by the proteasome to produce p50 and p52, respectively. Rel proteins bind p50 and p52 to form dimeric complexes that bind DNA and regulate transcription. In unstimulated cells, NF-kB is sequestered in the cytoplasm by IκB inhibitory proteins. NF-κB-activating agents can induce the phosphorylation of IkB proteins, targeting them for rapid degradation through the ubiquitin-proteasome pathway and releasing NF-kB to enter the nucleus where it regulates gene expression. NIK and IKKα (IKK1) regulate the phosphorylation and processing of NF-κB2 (p100) to produce p52, which translocates to the nucleus.

#### **Product:**

Liquid in 0.42% Potassium phosphate, 0.87% Sodium chloride, pH 7.3, 30% glycerol, and 0.01% sodium azide.

#### **Molecular Weight:**

~ 65 kDa

#### **Swiss-Prot:**

Q04206

#### **Purification&Purity:**

The antibody was affinity-purified from mouse antiserum by affinity-chromatography using epitope-specific immunogen and the purity is > 95% (by SDS-PAGE).

## **Applications:**

WB (1/500 - 1/1000), IHC (1/200 - 1/500), IF/ICC (1/100 - 1/200), IP (1/100 - 1/200)

#### Storage&Stability:

Store at  $4\,\mathrm{C}$  short term. Aliquot and store at -20  $\mathrm{C}$  long term. Avoid freeze-thaw cycles.

### **Specificity:**

## Bioworld Technology, Inc.


Add: 1660 South Highway 100, Suite 500 St. Louis Park,

MN 55416,USA.


Email: <u>info@bioworlde.com</u>

Tel: 6123263284 Fax: 6122933841 Recognizes endogenous levels of NF-kappaB p65 protein.


#### **DATA:**

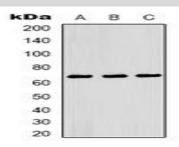

Western blot analysis of NF-kappaB p65 expression in Hela (A), mouse spleen (B), rat heart (C) whole cell lysates.

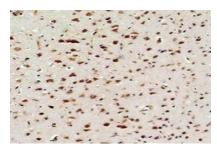

Immunohistochemical analysis of NF-kappaB p65 staining in rat brain formalin fixed paraffin embedded tissue section. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH 6.0). The section was then incubated with the antibody at room temperature and detected using an HRP conjugated compact polymer system. DAB was used as the chromogen. The section was then counterstained with haematoxylin and mounted with DPX.

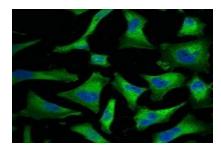

Immunofluorescent analysis of NF-kappaB p65 staining in Hela cells. Formalin-fixed cells were permeabilized with 0.1% Triton X-100 in TBS for 5-10 minutes and blocked with 3% BSA-PBS for 30 minutes at

## Bioworld technology, co. Ltd.

Add: No 9, weidi road Qixia District Nanjing, 210046,

P. R. China.


Email: <a href="mailto:info@biogot.com">info@biogot.com</a>
Tel: 0086-025-68037686
Fax: 0086-025-68035151



## PRODUCT DATA SHEET

Bioworld Technology,Inc.

room temperature. Cells were probed with the primary antibody in 3% BSA-PBS and incubated overnight at 4  $\,^{\circ}$ C in a hidified chamber. Cells were washed with PBST and incubated with a FITC-conjugated secondary antibody (green) in PBS at room temperature in the dark. DAPI was used to stain the cell nuclei (blue).

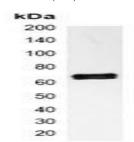

Immunoprecipitation of NF-kappaB p65 from 0.5mg Hela whole cell extract lysate, using Anti-NF-kappaB p65 Antibody.

#### Note:

For research use only, not for use in diagnostic procedure.

Bioworld Technology, Inc.

Add: 1660 South Highway 100, Suite 500 St. Louis Park,

MN 55416,USA.

Email: info@bioworlde.com

Tel: 6123263284 Fax: 6122933841 Bioworld technology, co. Ltd.

Add: No 9, weidi road Qixia District Nanjing, 210046,

P. R. China.

Email: <u>info@biogot.com</u>
Tel: 0086-025-68037686
Fax: 0086-025-68035151