

CD3e monoclonal antibody

Catalog: MB66179

Host:

Mouse

Reactivity: Human

BackGround:

When T cells encounter antigens via the T cell receptor (TCR), information about the quantity and quality of antigens is relayed to the intracellular signal transduction machinery. This activation process depends mainly on CD3 (Cluster of Differentiation 3), a multiunit protein complex that directly associates with the TCR. CD3 is composed of four polypeptides: ζ , γ , ϵ , and δ . Each of these polypeptides contains at least one immunoreceptor tyrosine-based activation motif (ITAM). Engagement of the TCR complex with foreign antigens induces tyrosine phosphorylation in the ITAM motifs and phosphorylated ITAMs function as docking sites for signaling molecules such as ZAP-70 and the p85 subunit of PI-3 kinase . TCR ligation also induces a conformational change in CD3_ε, such that a proline region is exposed and then associates with the adaptor protein Nck .

Product:

Mouse IgG2a. Liquid in PBS containing 50% glycerol, 0.2% BSA and 0.01% sodium azide.

Molecular Weight:

~ 24 kDa

Swiss-Prot:

P07766

Purification&Purity:

The antibody was affinity-purified from mouse antiserum by affinity-chromatography using epitope-specific immunogen and the purity is > 95% (by SDS-PAGE).

Applications:

WB (1/500 - 1/1000), IHC (1/100 - 1/300)

Storage&Stability:

Store at $4 \,^{\circ}{\rm C}$ short term. Aliquot and store at $-20 \,^{\circ}{\rm C}$ long term. Avoid freeze-thaw cycles.


Specificity:

Recognizes endogenous levels of CD3e protein.

DATA:



Immunohistochemical analysis of CD3e staining in human diffuse large b-cell lymphoma formalin fixed paraffin embedded tissue section. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH 6.0). The section was then incubated with the antibody at room temperature and detected using an HRP conjugated compact polymer system. DAB was used as the chromogen. The section was then counterstained with haematoxylin and mounted with DPX.

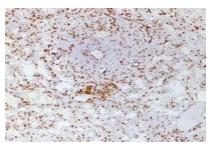

Immunohistochemical analysis of CD3e staining in human spleen formalin fixed paraffin embedded tissue section. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH 6.0). The section was then incubated with the antibody at room temperature and detected using an HRP conjugated compact polymer system. DAB was used as the chromogen. The section was then counterstained with haematoxylin and mounted with DPX.

Bioworld Technology, Inc.


Add:	1660 South Highway 100, Suite 500 St. Louis Park,
	MN 55416,USA.
Email:	info@bioworlde.com
Tel:	6123263284
Fax:	6122933841

Bioworld technology, co. Ltd. Add: No 9, weidi road Qixia District Nanjing, 210046, P. R. China. Email: <u>info@biogot.com</u> Tel: 0086-025-68037686 Fax: 0086-025-68035151



PRODUCT DATA SHEET

Bioworld Technology, Inc.

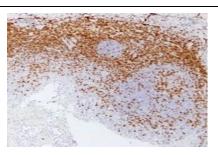

Immunohistochemical analysis of CD3e staining in human tonsil forma-

lin fixed paraffin embedded tissue section. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH 6.0). The section was then incubated with the antibody at room temperature and detected using an HRP conjugated compact polymer system. DAB was used as the chromogen. The section was then counterstained with haematoxylin and mounted with DPX.

Note:

For research use only, not for use in diagnostic procedure.

Bioworld Technology, Inc.		Bioworld technology, co. Ltd.	
Add:	1660 South Highway 100, Suite 500 St. Louis Park,	Add:	No 9, weidi road Qixia District Nanjing, 210046,
	MN 55416,USA.		P. R. China.
Email:	<u>info@bioworlde.com</u>	Email:	info@biogot.com
Tel:	6123263284	Tel:	0086-025-68037686
Fax:	6122933841	Fax:	0086-025-68035151