

## PRODUCT DATA SHEET

Bioworld Technology,Inc.

# Focal Adhesion Kinase Rabbit monoclonal antibody

Catalog: MB66341 Host: Rabbit Reactivity: Human, Mouse, Rat

#### **BackGround:**

Focal adhesion kinase (FAK) is a widely expressed cytoplasmic protein tyrosine kinase involved in integrin-mediated signal transduction. It plays an important role in the control of several biological processes, including cell spreading, migration, and survival. Activation of FAK by integrin clustering leads to autophosphorylation at Tyr397, which is a binding site for the Src family kinases PI3K and PLCγ. Recruitment of Src family kinases results in the phosphorylation of Tyr407, Tyr576, and Tyr577 in the catalytic domain, and Tyr871 and Tyr925 in the carboxy-terminal region of FAK.

## **Product:**

Liquid in 50mM Tris-Glycine (pH 7.4), 0.15M NaCl, 50% Glycerol, 0.01% Sodium azide and 0.05% BSA.

#### **Molecular Weight:**

~ 125 kDa

#### **Swiss-Prot:**

O05397

## **Purification&Purity:**

The antibody was purified by immunogen affinity chromatography.

#### **Applications:**

WB (1/500 - 1/1000), IHC (1/50 - 1/100), IF/ICC (1/50 - 1/100), IP (1/10 - 1/50)


#### **Storage&Stability:**

Store at  $4\,\mathrm{C}$  short term. Aliquot and store at  $-20\,\mathrm{C}$  long term. Avoid freeze-thaw cycles.


## **Specificity:**

Recognizes endogenous levels of Focal Adhesion Kinase protein.


#### **DATA:**

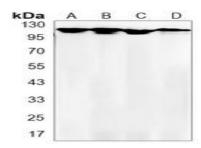

Western blot analysis of Focal Adhesion Kinase expression in K562 (A), C6 (B), NIH3T3 (C), Hela (D) whole cell lysates.

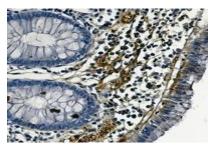

Immunohistochemical analysis of Focal Adhesion Kinase staining in human colon cancer formalin fixed paraffin embedded tissue section. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH 6.56). The section was then incubated with the antibody at room temperature and detected using an HRP conjugated compact polymer system. DAB was used as the chromogen. The section was then counterstained with haematoxylin and mounted with DPX.

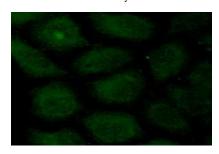

Immunofluorescent analysis of Focal Adhesion Kinase staining in A549 cells. Formalin-fixed cells were permeabilized with 0.1% Triton X-100 in TBS for 5-10 minutes and blocked with 3% BSA-PBS for 30 minutes at room temperature. Cells were probed with the primary antibody in 3% BSA-PBS and incubated overnight at 4 °C in a hidified chamber. Cells were washed with PBST and incubated with a AF488-conjugated

## Bioworld Technology, Inc.

Add: 1660 South Highway 100, Suite 500 St. Louis Park,

MN 55416,USA.

Email: info@bioworlde.com

Tel: 6123263284 Fax: 6122933841

## Bioworld technology, co. Ltd.

Add: No 9, weidi road Qixia District Nanjing, 210046,

P. R. China.

Email: <u>info@biogot.com</u> Tel: 0086-025-68037686 Fax: 0086-025-68035151



## PRODUCT DATA SHEET

Bioworld Technology,Inc.

secondary antibody (green) in PBS at room temperature in the dark.

For research use only, not for use in diagnostic procedure.

Note:

Bioworld Technology, Inc.

Add: 1660 South Highway 100, Suite 500 St. Louis Park,

MN 55416,USA.

Email: info@bioworlde.com

Tel: 6123263284 Fax: 6122933841 Bioworld technology, co. Ltd.

Add: No 9, weidi road Qixia District Nanjing, 210046,

P. R. China.

Email: <a href="mailto:info@biogot.com">info@biogot.com</a>
Tel: 0086-025-68037686
Fax: 0086-025-68035151