



# Catalase Microplate Assay Kit

## User Manual

Catalog # ASK1061

Detection and Quantification of Catalase (CAT) Activity in Urine, Serum, Plasma, Tissue extracts, Cell lysate, Cell culture media and Other biological fluids Samples.

**For research use only. Not for diagnostic or therapeutic procedures.**

Bioworld Technology, Inc. (USA)

Email: [info@bioworld.com](mailto:info@bioworld.com)

Web: [www.bioworld.com](http://www.bioworld.com)

Bioworld technology, co. Ltd. (China)

Email: [info@biogot.com](mailto:info@biogot.com)

Web: [www.biogot.com](http://www.biogot.com)



I. INTRODUCTION.....2

II. KIT COMPONENTS.....3

III. MATERIALS REQUIRED BUT NOT PROVIDED.....3

IV. SAMPLE PREPARATION.....4

V. ASSAY PROCEDURE.....4

VI. CALCULATION.....5

VII. TYPICAL DATA.....6



## **I. INTRODUCTION**

Catalase is an antioxidant enzyme ubiquitously present in mammalian and non-mammalian aerobic cells containing a cytochrome system. It was initially isolated from ox liver and later from blood, bacterial, and plant sources.

The enzyme contains 4 ferrihemoprotein groups per molecule. The enzyme has a molecular mass of 240 kDa. Catalase activity varies greatly between tissues. The activity is highest in the liver and kidney, and lowest in connective tissues. In eukaryotic cells the enzyme is concentrated in the subcellular organelles called peroxisomes (microbodies). Catalase catalyses the decomposition of hydrogen peroxide (H<sub>2</sub>O<sub>2</sub>) to water and oxygen. Hydrogen peroxide is formed in the eukaryotic cell as a by-product of various oxidase and superoxide dismutase reactions. Hydrogen peroxide is highly deleterious to the cell and its accumulation causes oxidation of cellular targets such as DNA, proteins, and lipids leading to mutagenesis and cell death. Removal of the H<sub>2</sub>O<sub>2</sub> from the cell by catalase provides protection against oxidative damage to the cell. It's role in oxidative stress related diseases has been widely studied.

The assay is initiated with the enzymatic hydrolysis of H<sub>2</sub>O<sub>2</sub> by CAT. The reaction product can be react with the dry reagent, and measured at a colorimetric readout at 405 nm.

**II. KIT COMPONENTS**

Component	Volume	Storage
96-Well Microplate	1 plate	
Assay Buffer	30 ml x 4	4 °C
Substrate Diluent	8 ml x 1	4 °C
Substrate	75 µl x 1	4 °C, keep in dark
Dye Reagent	Powder x 1	4 °C
Standard (100 mmol/L)	1 ml x 1	4 °C
Positive Control	Powder x 1	-20 °C
Technical Manual	1 Manual	

**Note:**

**Substrate:** add all Substrate into Substrate Diluent, mix.

**Dye Reagent:** add 10 ml distilled water to dissolve before use.

**Positive Control:** add 1 ml distilled water to dissolve before use.

**III. MATERIALS REQUIRED BUT NOT PROVIDED**

1. Microplate reader to read absorbance at 405 nm
2. Distilled water
3. Pipettor
4. Pipette tips
5. Mortar
6. Centrifuge
7. Timer
8. Ice



**IV. SAMPLE PREPARATION**

1. For cell and bacteria samples

Collect cell or bacteria into centrifuge tube, discard the supernatant after centrifugation, add 1 ml Assay buffer for  $5 \times 10^6$  cell or bacteria, sonicate (with power 20%, sonication 3s, interval 10s, repeat 30 times); centrifuged at 8000g 4 °C for 20 minutes, take the supernatant into a new centrifuge tube and keep it on ice for detection.

2. For tissue samples

Weigh out 0.1 g tissue, homogenize with 1 ml Assay buffer on ice, centrifuged at 8000g 4 °C for 20 minutes, take the supernatant into a new centrifuge tube and keep it on ice for detection.

3. For serum or plasma samples

Detect directly.

**V. ASSAY PROCEDURE**

Warm the Substrate to room temperature before use.

Add following reagents into the microplate:

Reagent	Sample	Blank	Standard	Positive Control
Substrate	80 $\mu$ l	80 $\mu$ l	--	80 $\mu$ l
Sample	20 $\mu$ l	--	--	--
Distilled water	--	20 $\mu$ l		--
Positive Control	--	--	--	20 $\mu$ l
Incubate at room temperature for 3 minutes.				
Standard	--	--	100 $\mu$ l	--



## PRODUCT DATA SHEET

Bioworld Technology, Inc

Dye Reagent	100 $\mu$ l	100 $\mu$ l	100 $\mu$ l	100 $\mu$ l
Mix, measured at 405 nm and record the absorbance.				

### VI. CALCULATION

**Unit Definition:** One unit of catalase activity is defined as the enzyme decomposes 1  $\mu$ mol of hydrogen peroxide per minute.

1. According to the protein concentration of sample

$$\begin{aligned} \text{CAT (U/mg)} &= (C_{\text{Standard}} \times V_{\text{Standard}}) \times (\text{OD}_{\text{Blank}} - \text{OD}_{\text{Sample}}) / \text{OD}_{\text{Standard}} / (V_{\text{Sample}} \times C_{\text{Protein}}) / T \\ &= 166.7 \times (\text{OD}_{\text{Blank}} - \text{OD}_{\text{Sample}}) / \text{OD}_{\text{Standard}} / C_{\text{Protein}} \end{aligned}$$

2. According to the weight of sample

$$\begin{aligned} \text{CAT (U/g)} &= (C_{\text{Standard}} \times V_{\text{Standard}}) \times (\text{OD}_{\text{Blank}} - \text{OD}_{\text{Sample}}) / \text{OD}_{\text{Standard}} / (W \times V_{\text{Sample}} / V_{\text{Assay}}) / T \\ &= 166.7 \times (\text{OD}_{\text{Blank}} - \text{OD}_{\text{Sample}}) / \text{OD}_{\text{Standard}} / W \end{aligned}$$

3. According to the quantity of cell or bacteria

$$\begin{aligned} \text{CAT (U}/10^4) &= (C_{\text{Standard}} \times V_{\text{Standard}}) \times (\text{OD}_{\text{Blank}} - \text{OD}_{\text{Sample}}) / \text{OD}_{\text{Standard}} / (N \times V_{\text{Sample}} / V_{\text{Assay}}) / T \\ &= 166.7 \times (\text{OD}_{\text{Blank}} - \text{OD}_{\text{Sample}}) / \text{OD}_{\text{Standard}} / N \end{aligned}$$

4. According to the volume of sample

$$\begin{aligned} \text{CAT (U/ml)} &= (C_{\text{Standard}} \times V_{\text{Standard}}) \times (\text{OD}_{\text{Blank}} - \text{OD}_{\text{Sample}}) / \text{OD}_{\text{Standard}} / V_{\text{Sample}} / T \\ &= 166.7 \times (\text{OD}_{\text{Blank}} - \text{OD}_{\text{Sample}}) / \text{OD}_{\text{Standard}} \end{aligned}$$

$C_{\text{Protein}}$ : the protein concentration, mg/ml;

$C_{\text{Standard}}$ : the concentration of Standard, 100 mmol/L = 100  $\mu$ mol/ml;

W: the weight of sample, g;

N: the quantity of cell or bacteria,  $N \times 10^4$ ;

$V_{\text{Sample}}$ : the volume of sample, 0.02 ml;


$V_{\text{Standard}}$ : the volume of sample, 0.1 ml;

$V_{\text{Assay}}$ : the volume of Assay buffer, 1 ml;


T: the reaction time, 3 minutes.

## VII. TYPICAL DATA

The standard curve is for demonstration only. A standard curve must be run with each assay.

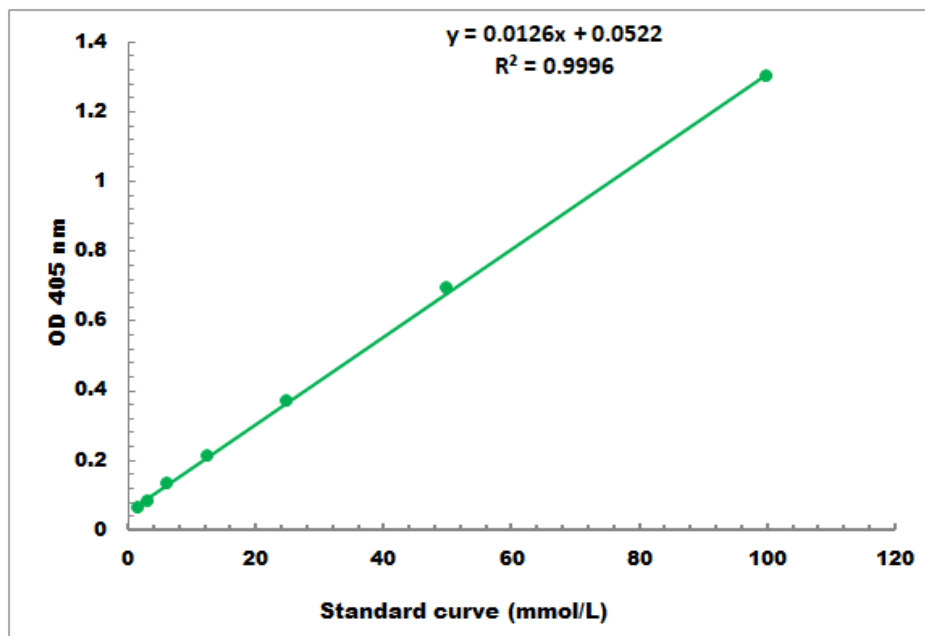

Detection Range: 1 mmol/L - 100 mmol/L

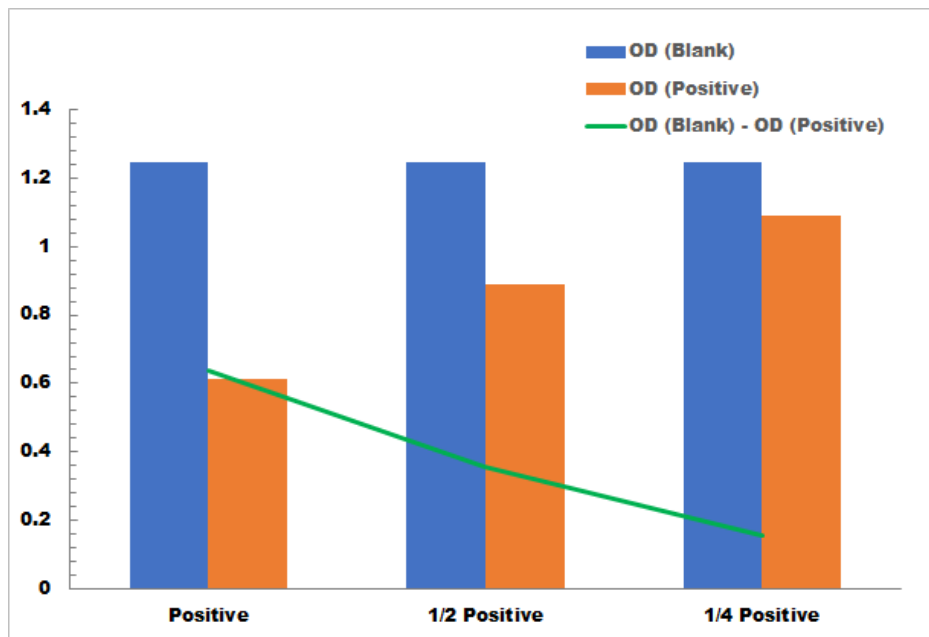

Positive Control reaction in 96-well plate assay with decreasing Positive Control concentration