



Kalium Microplate Assay Kit

User Manual

Catalog # ASK1104

Detection and Quantification of Kalium (K^+) Content in Serum, Urine, Saliva and Other biological fluids Samples.

For research use only. Not for diagnostic or therapeutic procedures.

Bioworld Technology, Inc. (USA)

Email: info@bioworldde.com

Web: www.bioworldde.com

Bioworld technology, co. Ltd. (China)

Email: info@biogot.com

Web: www.biogot.com



I. INTRODUCTION.....	2
II. KIT COMPONENTS.....	3
III. MATERIALS REQUIRED BUT NOT PROVIDED.....	3
IV. SAMPLE PREPARATION.....	4
V. ASSAY PROCEDURE.....	5
VI. CALCULATION.....	6
VII. TYPICAL DATA.....	7



I. INTRODUCTION

Kalium is necessary for the function of all living cells. Kalium shifts across nerve cell membranes are necessary for normal nerve transmission: kalium depletion can result in numerous abnormalities, including of heart rhythm. Fresh fruits and vegetables are good dietary sources of kalium.

The kalium can react with sodium tetraphenylborate. The products kalium tetraphenylboron can be measured at a colorimetric readout at 520 nm.

**II. KIT COMPONENTS**

Component	Volume	Storage
96-Well Microplate	1 plate	
Assay Buffer	25 ml x 2	4 °C
Reaction Buffer	4 ml x 1	4 °C
Dye Reagent	12 ml x 1	4 °C
Standard (2 mmol/L)	1 ml x 1	4 °C
Technical Manual	1 Manual	

III. MATERIALS REQUIRED BUT NOT PROVIDED

1. Microplate reader to read absorbance at 520 nm
2. Distilled water
3. Pipettor
4. Pipette tips
5. Mortar
6. Centrifuge
7. Timer



IV. SAMPLE PREPARATION

1. For serum sample

Add 50 μ l serum and 450 μ l Assay buffer into the microcentrifuge tube, mix, centrifuged at 8,000g 25 °C for 10 minutes, take the supernatant into a new centrifuge tube for detection.

**V. ASSAY PROCEDURE**

Add following reagents into the microplate:

Reagent	Sample	Standard	Blank
Sample	40 μ l	--	--
Standard	--	40 μ l	--
Distilled water	--	--	40 μ l
Reaction Buffer	40 μ l	40 μ l	40 μ l
Mix, wait for 5 minutes.			
Dye Reagent	120 μ l	120 μ l	120 μ l
Mix, measured at 520 nm and record the absorbance.			



VI. CALCULATION


1. According to the serum sample

$$\begin{aligned} K^+ (\text{mmol/L}) &= C_{\text{Standard}} \times (\text{OD}_{\text{Sample}} - \text{OD}_{\text{Blank}}) / (\text{OD}_{\text{Standard}} - \text{OD}_{\text{Blank}}) \times 10 \\ &= 20 \times (\text{OD}_{\text{Sample}} - \text{OD}_{\text{Blank}}) / (\text{OD}_{\text{Standard}} - \text{OD}_{\text{Blank}}) \end{aligned}$$

C_{Standard}: the concentration of Standard, 2 mmol/L.

**VII. TYPICAL DATA**

The standard curve is for demonstration only. A standard curve must be run with each assay.

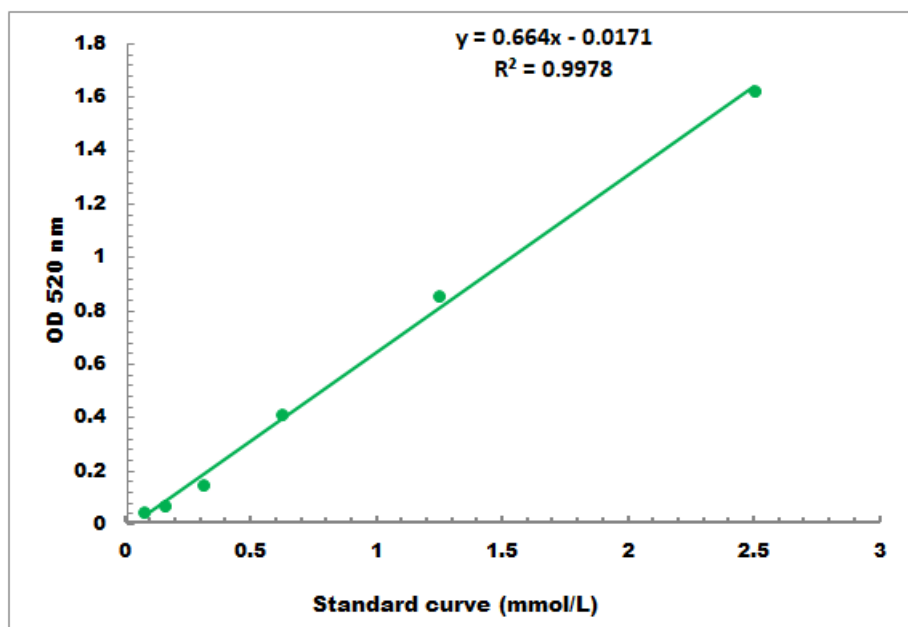

Detection Range: 0.2 mmol/L - 2 mmol/L