



Zinc Microplate Assay Kit

User Manual

Catalog # ASK1110

Detection and Quantification of Zinc (Zn^{2+}) Content in Serum, Urine, Saliva and Other biological fluids Samples.

For research use only. Not for diagnostic or therapeutic procedures.

Bioworld Technology, Inc. (USA)

Email: info@bioworldde.com

Web: www.bioworldde.com

Bioworld technology, co. Ltd. (China)

Email: info@biogot.com

Web: www.biogot.com



I. INTRODUCTION.....2

II. KIT COMPONENTS.....3

III. MATERIALS REQUIRED BUT NOT PROVIDED.....3

IV. SAMPLE PREPARATION.....4

V. ASSAY PROCEDURE.....5

VI. CALCULATION.....6

VII. TYPICAL DATA.....7



I. INTRODUCTION

Zinc, a metallic chemical element, symbol Zn and atomic number 30 is chemically similar to Magnesium due to its similar size and sole oxidation state of 2^+ . Zinc is an essential mineral of great biological significance, because many enzymes require it as an essential cofactor. Examples of zinc's biological roles include signal transduction, gene expression, regulation of apoptosis, synaptic plasticity and prostate gland function.

The reaction products can be measured at a colorimetric readout at 558 nm.



II. KIT COMPONENTS

Component	Volume	Storage
96-Well Microplate	1 plate	
Assay Buffer	30 ml x 1	4 °C
Dye Reagent	Powder x 1	4 °C
Reaction Buffer	20 ml x 1	4 °C
Standard (100 µmol/L)	1 ml x 1	4 °C
Technical Manual	1 Manual	

Note:

Dye Reagent: add 0.2 ml alcohol to dilute before use.

Working solution: add all Dye Reagent into the Reaction Buffer, mix.

III. MATERIALS REQUIRED BUT NOT PROVIDED

1. Microplate reader to read absorbance at 558 nm
2. Distilled water
3. Pipettor
4. Pipette tips
5. Alcohol
6. Centrifuge
7. Timer



IV. SAMPLE PREPARATION

1. For serum sample

Add 100 μ l serum and 100 μ l assay buffer into the microcentrifuge tube, mix, centrifuged at 10,000g 4 °C for 10 minutes, take the supernatant into a new centrifuge tube for detection.



V. ASSAY PROCEDURE

Add following reagents into the microplate:

Reagent	Blank	Standard	Sample
Distilled water	25 μ l	--	--
Assay Buffer	25 μ l	25 μ l	--
Standard	--	25 μ l	--
Sample	--	--	50 μ l
Working solution	200 μ l	200 μ l	200 μ l
Mix, wait for 2 minutes, measured at 558 nm and record the absorbance.			



VI. CALCULATION


1. According to the serum sample

$$\begin{aligned} \text{Zn}^{2+} (\mu\text{mol/L}) &= C_{\text{Standard}} \times (\text{OD}_{\text{Sample}} - \text{OD}_{\text{Blank}}) / (\text{OD}_{\text{Standard}} - \text{OD}_{\text{Blank}}) \\ &= 100 \times (\text{OD}_{\text{Sample}} - \text{OD}_{\text{Blank}}) / (\text{OD}_{\text{Standard}} - \text{OD}_{\text{Blank}}) \end{aligned}$$

C_{Standard} : the concentration of Standard, 100 $\mu\text{mol/L}$.

VII. TYPICAL DATA

The standard curve is for demonstration only. A standard curve must be run with each assay.

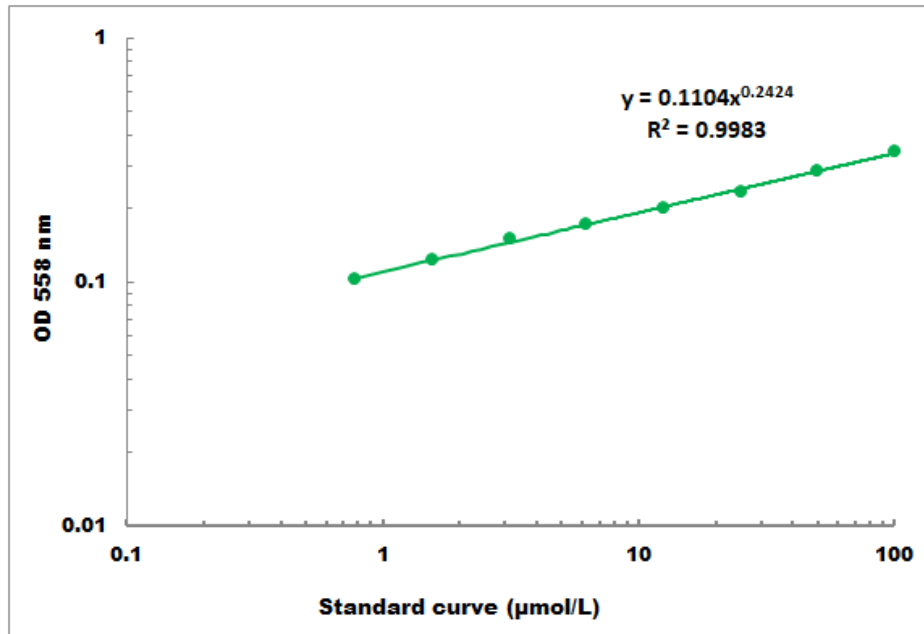

Detection Range: 1 µmol/L - 200 µmol/L