

PRODUCT DATA SHEET

Bioworld Technology,Inc.

Prohibitin polyclonal antibody

Catalog: BS6163 Host: Rabbit Reactivity: Human, Mouse, Rat

BackGround:

This gene is evolutionarily conserved, and its product is proposed to play a role in human cellular senescence and tumor suppression. Antiproliferative activity is reported to be localized to the 3' UTR, which is proposed to function as a trans-acting regulatory RNA. Several pseudogenes of this gene have been identified. Alternative splicing results in multiple transcript variants.

Product:

1mg/ml in PBS with 0.02% sodium azide, 50% glycerol, pH7.2

Molecular Weight:

30KDa

Swiss-Prot:

P35232

Purification&Purity:

The antibody was affinity-purified from rabbit antiserum by affinity-chromatography using epitope-specific immunogen and the purity is > 95% (by SDS-PAGE).

Applications:

WB,1:500 - 1:2000|IHC,1:50 - 1:200|IF/ICC,1:50 - 1:200


Storage&Stability:

Store at $4 \, \mathbb{C}$ short term. Aliquot and store at $-20 \, \mathbb{C}$ long term. Avoid freeze-thaw cycles.

Category:

Polyclonal Antibodies

DATA:

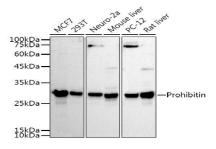

Western blot analysis of extracts of various cell lines, using Prohibitin

Rabbit pAb antibody at 1:500 dilution.


Secondary antibody: HRP Goat Anti-Rabbit IgG at 1:10000 dilution.

br/>Lysates/proteins: 25ug per lane.


br/>Blocking buffer: 3% nonfat dry milk in TBST.

br/>Detection: ECL Basic Kit .


cyp/Exposure time: 10s.

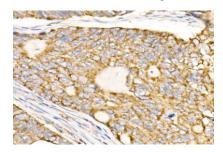

Immunohistochemistry of paraffin-embedded human colon carcinoma using Prohibitin Rabbit pAb at dilution of 1:25 .Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

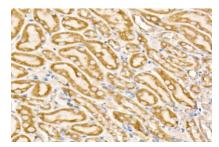

Immunohistochemistry of paraffin-embedded mouse kidney using Prohibitin Rabbit pAb at dilution of 1:25 .Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

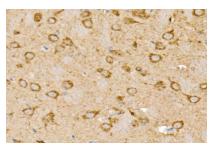

Immunohistochemistry of paraffin-embedded rat brain using Prohibitin Rabbit pAb at dilution of 1:25 .Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining

Bioworld Technology, Inc.

Add: 1660 South Highway 100, Suite 500 St. Louis Park,

MN 55416,USA.

Email: info@bioworlde.com

Tel: 6123263284 Fax: 6122933841

Bioworld technology, co. Ltd.

Add: No 9, weidi road Qixia District Nanjing, 210046,

P. R. China.

Email: info@biogot.com
Tel: 0086-025-68037686
Fax: 0086-025-68035151



PRODUCT DATA SHEET

Bioworld Technology,Inc.

neo	tocol	1

For research use only, not for use in diagnostic procedure.

Note:

Bioworld Technology, Inc.

Add: 1660 South Highway 100, Suite 500 St. Louis Park,

MN 55416,USA.

Email: info@bioworlde.com

Tel: 6123263284 Fax: 6122933841 Bioworld technology, co. Ltd.

Add: No 9, weidi road Qixia District Nanjing, 210046,

P. R. China.

Email: info@biogot.com
Tel: 0086-025-68037686
Fax: 0086-025-68035151