

PRODUCT DATA SHEET

Bioworld Technology,Inc.

Cyclin K polyclonal antibody

Catalog: BS72296 Host: Rabbit Reactivity: Human, Mouse, Rat

BackGround:

The protein encoded by this gene is a member of the transcription cyclin family. These cyclins may regulate transcription through their association with and activation of cyclin-dependent kinases (CDK) that phosphorylate the C-terminal domain (CTD) of the large subunit of RNA polymerase II. This gene product may play a dual role in regulating CDK and RNA polymerase II activities.

Product:

1 mg/ml in PBS with 0.02% sodium azide, 50% glycerol, pH7.2

Molecular Weight:

90kDa

Swiss-Prot:

O75909

Purification&Purity:

The antibody was affinity-purified from rabbit antiserum by affinity-chromatography using epitope-specific immunogen and the purity is > 95% (by SDS-PAGE).

Applications:

WB,1:500 - 1:2000|IHC,1:50 - 1:200


Storage&Stability:

Store at $4 \,\mathrm{C}$ short term. Aliquot and store at $-20 \,\mathrm{C}$ long term. Avoid freeze-thaw cycles.

Category:

Polyclonal Antibodies

DATA:

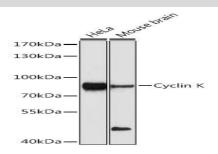

Western blot analysis of extracts of various cell lines, using Cyclin K

Antibody at 1:1000 dilution.


secondary antibody: HRP Goat Anti-Rabbit IgG at 1:10000 dilution.

br/>Lysates/proteins: 25ug per lane.


br/>Blocking buffer: 3% nonfat dry milk in TBST.

br/>Detection: ECL Enhanced Kit .


cypr/>Exposure time: 90s.

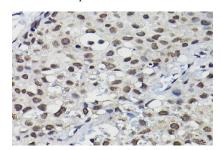

Immunohistochemistry of paraffin-embedded human breast cancer using Cyclin K Rabbit pAb at dilution of 1:100 .Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

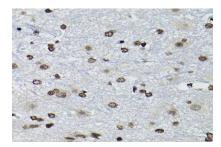

Immunohistochemistry of paraffin-embedded mouse spinal cord using Cyclin K Rabbit pAb at dilution of 1:100 .Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

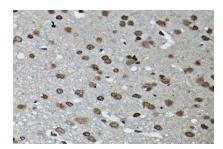

Immunohistochemistry of paraffin-embedded rat brain using Cyclin K Rabbit pAb at dilution of 1:100 .Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC

Bioworld Technology, Inc.

Add: 1660 South Highway 100, Suite 500 St. Louis Park,

MN 55416,USA.

Email: <u>info@bioworlde.com</u>

Tel: 6123263284 Fax: 6122933841

Bioworld technology, co. Ltd.

Add: No 9, weidi road Qixia District Nanjing, 210046,

P. R. China.

Email: <u>info@biogot.com</u>
Tel: 0086-025-68037686
Fax: 0086-025-68035151



PRODUCT DATA SHEET

Bioworld Technology,Inc.

staining protocol.

Note:

For research use only, not for use in diagnostic procedure.

Bioworld Technology, Inc.

Add: 1660 South Highway 100, Suite 500 St. Louis Park,

MN 55416,USA.

Email: <u>info@bioworlde.com</u>
Tel: 6123263284

Tel: 6123263284 Fax: 6122933841 Bioworld technology, co. Ltd.

Add: No 9, weidi road Qixia District Nanjing, 210046,

P. R. China.

Email: info@biogot.com
Tel: 0086-025-68037686
Fax: 0086-025-68035151