

GLB1 polyclonal antibody

Cata	log:	BS7306
------	------	--------

Host: 1

Rabbit

Reactivity: Human, Mouse, Rat

BackGround:

This gene encodes a member of the glycosyl hydrolase 35 family of proteins. Alternative splicing results in multiple transcript variants, at least one of which encodes a preproprotein that is proteolytically processed to generate the mature lysosomal enzyme. This enzyme catalyzes the hydrolysis of a terminal beta-linked galactose residue from ganglioside substrates and other glycoconjugates. Mutations in this gene may result in GM1-gangliosidosis and Morquio B syndrome.

Product:

1mg/ml in PBS with 0.02% sodium azide, 50% glycerol, pH7.2

Molecular Weight:

65KDa

Swiss-Prot:

P16278

Purification&Purity:

The antibody was affinity-purified from rabbit antiserum by affinity-chromatography using epitope-specific immunogen and the purity is > 95% (by SDS-PAGE).

Applications:

WB,1:500 - 1:2000|IHC,1:50 - 1:200|IF/ICC,1:50 - 1:200 Storage&Stability:

Store at $4 \,^{\circ}{\rm C}$ short term. Aliquot and store at $-20 \,^{\circ}{\rm C}$ long term. Avoid freeze-thaw cycles.


Category:

Polyclonal Antibodies


DATA:

Western blot analysis of extracts of various cell lines, using GLB1 antibody at 1:1000 dilution.
br/>Secondary antibody: HRP Goat Anti-Rabbit IgG at 1:10000 dilution.
br/>Lysates/proteins: 25ug per lane.
br/>Blocking buffer: 3% nonfat dry milk in TBST.
br/>Detection: ECL Basic Kit .< br/>Exposure time: 90s.

Immunohistochemistry of paraffin-embedded human kidney using GLB1 Rabbit pAb at dilution of 1:25 .Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

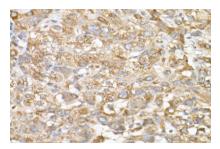

Immunohistochemistry of paraffin-embedded human liver cancer using GLB1 Rabbit pAb at dilution of 1:25 .Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

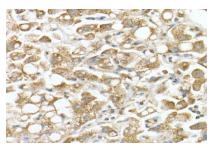

Immunohistochemistry of paraffin-embedded human liver using GLB1 Rabbit pAb at dilution of 1:25 .Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

Note:

For research use only, not for use in diagnostic procedure.

Bioworld Technology, Inc.

 Add:
 1660 South Highway 100, Suite 500 St. Louis Park, MN 55416,USA.

 Email:
 info@bioworlde.com

 Tel:
 6123263284

 Fax:
 6122933841

Bioworld technology, co. Ltd. Add: No 9, weidi road Qixia District Nanjing, 210046,

Add: No 9, weld road Qixia District Nanjing, 21004 P. R. China. Email: <u>info@biogot.com</u> Tel: 0086-025-68037686 Fax: 0086-025-68035151