

PRODUCT DATA SHEET

Bioworld Technology,Inc.

DBT polyclonal antibody

Catalog: BS77998 Host: Rabbit Reactivity: Human, Mouse, Rat

BackGround:

The branched-chain alpha-keto acid dehydrogenase complex (BCKD) is an inner-mitochondrial enzyme complex involved in the breakdown of the branched-chain amino acids isoleucine, leucine, and valine. The BCKD complex is thought to be composed of a core of 24 transacylase (E2) subunits, and associated decarboxylase (E1), dehydrogenase (E3), and regulatory subunits. This gene encodes the transacylase (E2) subunit. Mutations in this gene result in maple syrup urine disease, type 2. Alternatively spliced transcript variants have been described, but their biological validity has not been determined.

Product:

1 mg/ml in PBS with 0.02% sodium azide, 50% glycerol, pH7.2

Molecular Weight:

53KDa

Swiss-Prot:

P11182

Purification&Purity:

The antibody was affinity-purified from rabbit antiserum by affinity-chromatography using epitope-specific immunogen and the purity is > 95% (by SDS-PAGE).

Applications:

WB,1:500 - 1:2000|IHC,1:50 - 1:200|IF/ICC,1:50 - 1:200


Storage&Stability:

Store at $4 \,\mathrm{C}$ short term. Aliquot and store at $-20 \,\mathrm{C}$ long term. Avoid freeze-thaw cycles.

Modification:

Unmodification

DATA:

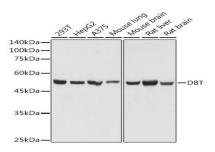

Western blot analysis of extracts of various cell lines, using DBT antibody at 1:1000 dilution.

Secondary antibody: HRP Goat Anti-Rabbit IgG at 1:10000 dilution.


br/>Lysates/proteins: 25ug per lane.

lane.

Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit.

Exposure time: 1s.

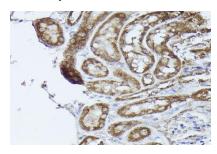

Immunohistochemistry of paraffin-embedded mouse kidney using DBT Rabbit pAb at dilution of 1:50 .Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

Immunofluorescence analysis of A-431 cells using DBT Rabbit pAb at dilution of 1:100 . Blue: DAPI for nuclear staining.

Immunofluorescence analysis of PC-12 cells using DBT Rabbit pAb

at dilution of 1:100 . Blue: DAPI for nuclear staining.

Note:

For research use only, not for use in diagnostic procedure.

Bioworld Technology, Inc.

Add: 1660 South Highway 100, Suite 500 St. Louis Park,

MN 55416,USA.

Email: <u>info@bioworlde.com</u>

Tel: 6123263284 Fax: 6122933841 Bioworld technology, co. Ltd.

Add: No 9, weidi road Qixia District Nanjing, 210046,

P. R. China.

Email: <u>info@biogot.com</u>
Tel: 0086-025-68037686
Fax: 0086-025-68035151