

PRODUCT DATA SHEET

Bioworld Technology,Inc.

MRPL13 polyclonal antibody

Catalog: BS78715 Host: Rabbit Reactivity: Human, Mouse, Rat

BackGround:

Mammalian mitochondrial ribosomal proteins are encoded by nuclear genes and help in protein synthesis within the mitochondrion. Mitochondrial ribosomes (mitoribosomes) consist of a small 28S subunit and a large 39S subunit. They have an estimated 75% protein to rRNA composition compared to prokaryotic ribosomes, where this ratio is reversed. Another difference between mammalian mitoribosomes and prokaryotic ribosomes is that the latter contain a 5S rRNA. Among different species, the proteins comprising the mitoribosome differ greatly in sequence, and sometimes in biochemical properties, which prevents easy recognition by sequence homology. This gene encodes a 39S subunit protein.

Product:

1 mg/ml in PBS with 0.02% sodium azide, 50% glycerol, pH7.2

Molecular Weight:

20kDa

Swiss-Prot:

Q9BYD1

Purification&Purity:

The antibody was affinity-purified from rabbit antiserum by affinity-chromatography using epitope-specific immunogen and the purity is > 95% (by SDS-PAGE).

Applications:

WB,1:500 - 1:2000|IHC,1:50 - 1:200


Storage&Stability:

Store at $4 \, \mathbb{C}$ short term. Aliquot and store at $-20 \, \mathbb{C}$ long term. Avoid freeze-thaw cycles.


Modification:

Unmodification


DATA:

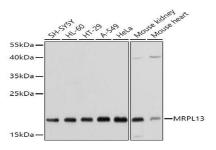

Western blot analysis of extracts of various cell lines, using MRPL13 antibody at 1:1000 dilution.
Secondary antibody: HRP Goat Anti-Rabbit IgG at 1:10000 dilution.
br/>Lysates/proteins: 25ug per lane.
br/>Blocking buffer: 3% nonfat dry milk in TBST.
Detection: ECL Basic Kit .
br/>Exposure time: 15s.

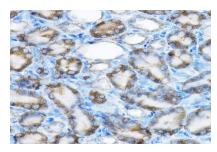

Immunohistochemistry of paraffin-embedded rat kidney using MRPL13 Rabbit pAb at dilution of 1:100 .Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

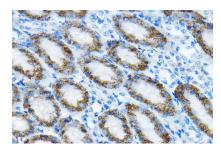

Immunohistochemistry of paraffin-embedded mouse kidney using MRPL13 Rabbit pAb at dilution of 1:100 .Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

Note:

For research use only, not for use in diagnostic procedure.

Bioworld Technology, Inc.

Add: 1660 South Highway 100, Suite 500 St. Louis Park,

MN 55416,USA.

Email: <u>info@bioworlde.com</u>

Tel: 6123263284 Fax: 6122933841

Bioworld technology, co. Ltd.

Add: No 9, weidi road Qixia District Nanjing, 210046,

P. R. China.

Email: <u>info@biogot.com</u>
Tel: 0086-025-68037686
Fax: 0086-025-68035151



PRODUCT DATA SHEET

Bioworld Technology,Inc.

Bioworld Technology, Inc.

Add: 1660 South Highway 100, Suite 500 St. Louis Park,

MN 55416,USA.

Email: info@bioworlde.com

Tel: 6123263284 Fax: 6122933841 Bioworld technology, co. Ltd.

Add: No 9, weidi road Qixia District Nanjing, 210046,

P. R. China.

Email: <u>info@biogot.com</u> Tel: 0086-025-68037686 Fax: 0086-025-68035151