

## PRODUCT DATA SHEET

Bioworld Technology,Inc.

# Ku70 polyclonal antibody

Catalog: BS78853 Host: Rabbit Reactivity: Human, Mouse, Rat

#### **BackGround:**

The p70/p80 autoantigen is a nuclear complex consisting of two subunits with molecular masses of approximately 70 and 80 kDa. The complex functions as a single-stranded DNA-dependent ATP-dependent helicase. The complex may be involved in the repair of nonhomologous DNA ends such as that required for double-strand break repair, transposition, and V(D)J recombination. High levels of autoantibodies to p70 and p80 have been found in some patients with systemic lupus erythematosus.

#### **Product:**

1mg/ml in PBS with 0.02% sodium azide, 50% glycerol, pH7.2

#### **Molecular Weight:**

70kDa

### **Swiss-Prot:**

P12956

#### **Purification&Purity:**

The antibody was affinity-purified from rabbit antiserum by affinity-chromatography using epitope-specific immunogen and the purity is > 95% (by SDS-PAGE).

#### **Applications:**

WB,1:500 - 1:2000|IHC,1:50 - 1:100|IF/ICC,1:50 - 1:100|IP,1:20 - 1:50


#### **Storage&Stability:**

Store at  $4\,\mathrm{C}$  short term. Aliquot and store at  $-20\,\mathrm{C}$  long term. Avoid freeze-thaw cycles.


#### **Modification:**

Unmodification


#### **DATA:**

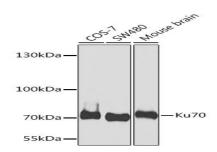

Western blot analysis of extracts of various cell lines, using Ku70 antibody at 1:1000 dilution.<br/>
Secondary antibody: HRP Goat Anti-Rabbit IgG at 1:10000 dilution.<br/>
br/>Lysates/proteins: 25ug per lane.<br/>
lane.<br/>
Blocking buffer: 3% nonfat dry milk in TBST.<br/>
Detection: ECL Basic Kit.<br/>
Exposure time: 60s.

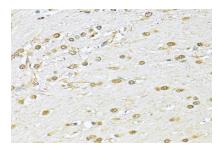

Immunohistochemistry of paraffin-embedded rat brain using Ku70 antibody at dilution of 1:100 .Perform microwave antigen retrieval with 10 mM PBS buffer pH 7.2 before commencing with IHC staining protocol.

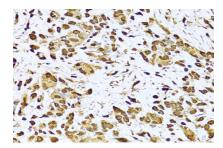

Immunohistochemistry of paraffin-embedded human gastric cancer using Ku70 antibody at dilution of 1:100 .Perform microwave antigen retrieval with 10 mM PBS buffer pH 7.2 before commencing with IHC staining protocol.

# Bioworld Technology, Inc.

Add: 1660 South Highway 100, Suite 500 St. Louis Park,

MN 55416,USA.

Email: info@bioworlde.com

Tel: 6123263284 Fax: 6122933841

# Bioworld technology, co. Ltd.

Add: No 9, weidi road Qixia District Nanjing, 210046,

P. R. China.

Email: <a href="mailto:info@biogot.com">info@biogot.com</a>
Tel: 0086-025-68037686
Fax: 0086-025-68035151



# PRODUCT DATA SHEET

Bioworld Technology,Inc.



Immunohistochemistry of paraffin-embedded human uterus using Ku70 antibody at dilution of 1:100 .Perform microwave antigen retrieval with 10 mM PBS buffer pH 7.2 before commencing with IHC staining protocol.

#### Note:

For research use only, not for use in diagnostic procedure.

# Bioworld Technology, Inc.

Add: 1660 South Highway 100, Suite 500 St. Louis Park,

MN 55416,USA.

Email: info@bioworlde.com

Tel: 6123263284 Fax: 6122933841 Bioworld technology, co. Ltd.

Add: No 9, weidi road Qixia District Nanjing, 210046,

P. R. China.

Email: <u>info@biogot.com</u>
Tel: 0086-025-68037686
Fax: 0086-025-68035151