



# Rat BDNF ELISA Kit User Manual

Catalog # CEK1577

Sandwich Enzyme-Linked Immunosorbent Assay for Quantitative

Detection of Rat BDNF Concentrations in Cell Culture Supernatants,

Serum, Plasma.

For research use only. Not for diagnostic or therapeutic procedures.

Bioworld Technology, Inc. (USA)

Email: <a href="mailto:info@bioworlde.com">info@bioworlde.com</a>
Web: <a href="mailto:www.bioworlde.com">www.bioworlde.com</a>

Bioworld technology, co. Ltd. (China)



I. INTRODUCTION	2
II. ASSAY PRINCIPLES	3
III. KIT COMPONENTS	4
IV. STORAGE AND STABILITY	4
V. MATERIALS REQUIRED BUT NOT PROVIDED	5
VI. HEALTH AND SAFETY PRECAUTIONS	5
VII. REAGENT PREPARATION	6
VIII. ASSAY PROCEDURE	9
IX. ASSAY PROCEDURE SUMMARY	11
X. TYPICAL DATA	12
XI. SENSITIVITY	12
XII. SPECIFICITY	12
XIII. CROSS REACTIVITY	13
XIV. REFERENCES	13
XV. TROUBLESHOOTING GUIDE	14



#### I. INTRODUCTION

Brain-derived neurotrophic factor (BDNF) is a prosurvival factor induced by cortical neurons that is necessary for survival of striatal neurons in the brain. It is a secreted protein with the molecular weight of 27.8kDa, consisting of 247 amino acids. It is known to promote neuronal survival and differentiation. BDNF shares substantial amino acid sequence identity with nerve growth factor (NGF). BDNF and neurotrophin-3 (NT-3) are two recently cloned neurotrophic factors that are homologous to NGF. mRNA products of the BDNF and NT-3 genes are detected in the adult human brain, suggesting that these proteins are involved in the maintenance of the adult nervous system. BDNF and other neurotrophins are critically involved in long-term potentiation (LTP). BDNF-mediated LTP is induced postsynaptically. BDNF has trophic effects on serotonergic (5-HT) neurons in the central nervous system. BDNF has an essential maintenance function in the regulation of anxiety-related behavior and in food intake through central mediators in both the basal and fasted state. It plays a role in treating breathing disorders such as respiratory insufficiency after spinal injury. The mature form of BDNF is identical in all mammals examined, and the gene encoding human BDNF to chromosome 11, band p13.

2

Bioworld technology, co. Ltd. (China)

Email: <a href="mailto:info@bioworlde.com">info@bioworlde.com</a>
Web: <a href="mailto:www.bioworlde.com">www.bioworlde.com</a>

Bioworld Technology, Inc. (USA)



#### II. ASSAY PRINCIPLES

The Cohesion Bioscience Rat BDNF ELISA (Enzyme-Linked Immunosorbent Assay) kit is an in vitro enzyme-linked immunosorbent assay for the quantitative measurement of Rat BDNF in Cell Culture Supernatants, Serum, Plasma. This assay employs an antibody specific for Rat BDNF coated on a 96-well plate. Standards and samples are pipetted into the wells and BDNF present in a sample is bound to the wells by the immobilized antibody. The wells are washed and biotinylated anti-Rat BDNF antibody is added. After washing away unbound biotinylated antibody, HRP-conjugated streptavidin is pipetted to the wells. The wells are again washed, a TMB substrate solution is added to the wells and color develops in proportion to the amount of BDNF bound. The Stop Solution changes the color from blue to yellow, and the intensity of the color is measured at 450 nm.

3

Bioworld Technology, Inc. (USA) Email: <a href="mailto:info@bioworlde.com">info@bioworlde.com</a> Web: <a href="https://www.bioworlde.com">www.bioworlde.com</a>

Email: <a href="mailto:info@biogot.com">info@biogot.com</a>
Web: <a href="mailto:www.biogot.com">www.biogot.com</a>

Bioworld technology, co. Ltd. (China)



#### III. KIT COMPONENTS

Component	Volume
96-well Plate Coated With Anti-Rat BDNF Antibody	12 x 8 Strips
Rat BDNF Standard	10 ng x 2
Biotin-Labeled Detection Antibody (100X)	120 μΙ
Streptavidin-HRP (100X)	120 μΙ
Standard/Sample Diluent	30 ml
Detection Antibody Diluent	12 ml
Streptavidin-HRP Diluent	12 ml
Wash Buffer (20X)	30 ml
TMB Substrate Solution	12 ml
Stop Solution	12 ml
Plate Adhesive Strips	3 Strips
Technical Manual	1 Manual

## IV. STORAGE AND STABILITY

All kit components are stable at 2 to 8 °C. Standard (recombinant protein) should be stored at -20 °C or -80 °C (recommended at -80 °C) after reconstitution. Opened Microplate Wells or reagents may be store for up to 1 month at 2 to 8 °C. Return unused wells to the pouch containing desiccant pack, reseal along entire edge. Note: the kit can be used within one year if the whole kit is stored at -20 °C. Avoid repeated freeze-thaw cycles.

4

Bioworld Technology, Inc. (USA)
Email: <a href="mailto:info@bioworlde.com">info@bioworlde.com</a>
Web: <a href="mailto:www.bioworlde.com">www.bioworlde.com</a>

Bioworld technology, co. Ltd. (China)

## V. MATERIALS REQUIRED BUT NOT PROVIDED

- 1. Microplate reader capable of measuring absorbance at 450 nm.
- 2. Adjustable pipettes and pipette tips to deliver 2  $\mu$ l to 1 ml volumes.
- 3. Adjustable 1-25 ml pipettes for reagent preparation.
- 4. 100 ml and 1 liter graduated cylinders.
- 5. Absorbent paper.
- 6. Distilled or deionized water.
- 7. Computer and software for ELISA data analysis.
- 8. Tubes to prepare standard or sample dilutions.

#### VI. HEALTH AND SAFETY PRECAUTIONS

- 1. Reagents provided in this kit may be harmful if ingested, inhaled or absorbed through the skin. Please carefully review the MSDS for each reagent before conducting the experiment.
- 2. Stop Solution contains 2 N Sulfuric Acid ( $H_2SO_4$ ) and is an extremely corrosive agent. Please wear proper eye, hand and face protection when handling this material. When the experiment is finished, be sure to rinse the plate with copious amounts of running water to dilute the Stop Solution prior to disposing the plate.

5

Bioworld Technology, Inc. (USA)
Email: <a href="mailto:info@bioworlde.com">info@bioworlde.com</a>
Web: <a href="mailto:www.bioworlde.com">www.bioworlde.com</a>

Bioworld technology, co. Ltd. (China)

#### VII. REAGENT PREPARATION

#### 1. Sample Preparation

Store samples to be assayed within 24 hours at 2-8°C. For long-term storage, aliquot and freeze samples at -20°C. Avoid repeated freeze-thaw cycles.

**Cell culture supernates**: Remove particulates by centrifugation, assay immediately or aliquot and store samples at -20°C.

**Serum**: Allow the serum to clot in a serum separator tube (about 4 hours) at room temperature. Centrifuge at approximately 1000 X g for 15 minutes. Analyze the serum immediately or aliquot and store samples at -20°C.

**Plasma**: Collect plasma using heparin or EDTA as an anticoagulant. Centrifuge for 15 minutes at 1500 X g within 30 minutes of collection. Assay immediately or aliquot and store samples at -20°C.

**Cell Lysates:** Collect cells and rinse cells with PBS. Homogenize and lyse cells throughly in lysate solution. Centrifuge cell lysates at approximately 10000 X g for 5 minutes to remove debris. Aliquots of the cell lysates were removed and assayed.

**Bone Tissue:** Extract demineralized bone samples in 4 M Guanidine-HCl and protease inhibitors. Dissolve the final sample in 2 M Guanidine-HCl.

**Tissue Homogenates:** Rinse tissue with PBS to remove excess blood, chopped into 1-2 mm pieces, and homogenize with a tissue homogenizer in PBS or in lysate solution, lysate solution: tissue net weight = 10ml : 1g (i.e. Add 10ml lysate solution to 1g tissue). Centrifuge at approximately 5000 X g for 5 minutes. Assay immediately or aliquot and store homogenates at -20°C. Avoid repeated freeze-thaw cycles.

**Urine**: Urinary samples should be cleared by centrifugation and then can be used directly without dilution. Storage at -20°C.

## 2. Rat BDNF Standard Preparation

6

Bioworld Technology, Inc. (USA) Email: <a href="mailto:info@bioworlde.com">info@bioworlde.com</a> Web: <a href="https://www.bioworlde.com">www.bioworlde.com</a> Bioworld technology, co. Ltd. (China) Email: info@biogot.com

Web: www.biogot.com

Reconstitute the lyophilized Rat BDNF Standard by adding 1 ml of Standard/Sample Diluent to make the 10000 pg/ml standard stock solution. Allow solution to sit at room temperature for 5 minutes, then gently vortex to mix completely. Use within one hour of reconstituting. Two tubes of the standard (10 ng per tube) are included in each kit. Use one tube for each experiment.

Perform 2-fold serial dilutions of the top standards to make the standard curve within the range of this assay (31.2 pg/ml - 2000 pg/ml) as below. Standard/Sample Dilution Buffer serves as the zero standard (0 pg/ml).

Standard	Add	Into
2000 pg/ml	200 μl of the Standard (10000 pg/ml)	800 μl of the Standard/Sample Diluent
1000 pg/ml	500 μl of the Standard (2000 pg/ml)	500 μl of the Standard/Sample Diluent
500 pg/ml	500 μl of the Standard (1000 pg/ml)	500 μl of the Standard/Sample Diluent
250 pg/ml	500 μl of the Standard (500 pg/ml)	500 μl of the Standard/Sample Diluent
125 pg/ml	500 μl of the Standard (250 pg/ml)	500 μl of the Standard/Sample Diluent
62.5 pg/ml	500 μl of the Standard (125 pg/ml)	500 μl of the Standard/Sample Diluent
31.25 pg/ml	500 μl of the Standard (62.5 pg/ml)	500 μl of the Standard/Sample Diluent
0 ng/ml	1 ml of the Standard/Sample Diluent	

**Note:** The standard solutions are best used within 2 hours. The 10000 pg/ml standard solution should be stored at 4°C for up to 12 hours, or at -20°C for up to 48 hours. Avoid repeated freeze-thaw cycles.

3. Biotin-Labeled Detection Antibody Working Solution Preparation

The Biotin-Labeled Detection Antibody should be diluted in 1:100 with the Detection

Antibody Diluent and mixed thoroughly. The solution should be prepared no more than 2 hours prior to the experiment.

7

Bioworld Technology, Inc. (USA)
Email: <a href="mailto:info@bioworlde.com">info@bioworlde.com</a>
Web: <a href="mailto:www.bioworlde.com">www.bioworlde.com</a>

Bioworld technology, co. Ltd. (China) Email: info@biogot.com

Web: www.biogot.com



## 4. Streptavidin-HRP Working Solution Preparation

The Streptavidin-HRP should be diluted in 1:100 with the Streptavidin-HRP Diluent and mixed thoroughly. The solution should be prepared no more than 1 hour prior to the experiment.

## 5. Wash Buffer Working Solution Preparation

Pour entire contents (30 ml) of the Wash Buffer Concentrate into a clean 1,000 ml graduated cylinder. Bring final volume to 600 ml with glass-distilled or deionized water (1:20).

8

Bioworld Technology, Inc. (USA) Email: <u>info@bioworlde.com</u>

Web: www.bioworlde.com

Bioworld technology, co. Ltd. (China)



VIII. ASSAY PROCEDURE

The Streptavidin-HRP Working Solution and TMB Substrate Solution must be kept warm at 37°C for 30 minutes before use. When diluting samples and reagents, they must be mixed completely and evenly. Standard detection curve should be prepared for each experiment. The user will decide sample dilution fold by crude estimation of

protein amount in samples.

1. Add 100  $\mu$ l of each standard and sample into appropriate wells.

2. Cover well and incubate for 90 minutes at room temperature or over night at 4°C

with gentle shaking.

3. Remove the cover, discard the solution and wash plate 3 times with Wash Buffer

Working Solution, and each time let Wash Buffer Working Solution stay in the wells

for 1 - 2 minutes. Blot the plate onto paper towels or other absorbent material. Do

NOT let the wells completely dry at any time.

4. Add 100 μl of Biotin-Labeled Detection Antibody Working Solution into each well

and incubate the plate at 37°C for 60 minutes.

5. Wash plate 3 times with Wash Buffer Working Solution, and each time let Wash

Buffer Working Solution stay in the wells for 1 - 2 minutes. Discard the Wash Buffer

Working Solution and blot the plate onto paper towels or other absorbent material.

6. Add 100 µl of Streptavidin-HRP Working Solution into each well and incubate the

plate at 37°C for 45 minutes.

7. Wash plate 5 times with Wash Buffer Working Solution, and each time let wash

buffer stay in the wells for 1 - 2 minutes. Discard the wash buffer and blot the plate

onto paper towels or other absorbent material.

8. Add 100  $\mu$ l of TMB Substrate Solution into each well and incubate plate at 37°C in

dark for 30 minutes.

9

Bioworld Technology, Inc. (USA)

Email: info@bioworlde.com

Web: www.bioworlde.com

Bioworld technology, co. Ltd. (China)



bioworld PRODUCT DATA SHEET

9. Add 100  $\mu l$  of Stop Solution into each well. The color changes into yellow immediately.

10. Read the O.D. absorbance at 450nm in a microplate reader within 30 minutes after adding the Stop Solution.

For calculation, (the relative O.D.450) = (the O.D.450 of each well) - (the O.D.450 of Zero well). The standard curve can be plotted as the relative O.D.450 of each standard solution (Y) vs. the respective concentration of the standard solution (X). The concentration of the samples can be interpolated from the standard curve.

**Note:** If the samples measured were diluted, multiply the dilution factor to the concentrations from interpolation to obtain the concentration before dilution.

10

Bioworld Technology, Inc. (USA) Email: <a href="mailto:info@bioworlde.com">info@bioworlde.com</a> Web: <a href="https://www.bioworlde.com">www.bioworlde.com</a>

Email: <a href="mailto:info@biogot.com">info@biogot.com</a>
Web: <a href="mailto:www.biogot.com">www.biogot.com</a>

Bioworld technology, co. Ltd. (China)

#### IX. ASSAY PROCEDURE SUMMARY

- Prepare all reagents, samples and standards
- Add 100 μl Standard or Sample
- Wash plate 3 times with Wash Buffer Working Solution
- Add 100 µl Biotin-Labeled Detection Antibody Working Solution
- Wash plate 3 times with Wash Buffer Working Solution
- Add 100 μl Streptavidin-HRP Working Solution
- Wash plate 5 times with Wash Buffer Working Solution
- Add 100 µl TMB Substrate Solution
- Add 100 µl Stop Solution
- Read the plate at 450nm

11


Bioworld Technology, Inc. (USA)
Email: <a href="mailto:info@bioworlde.com">info@bioworlde.com</a>
Web: <a href="mailto:www.bioworlde.com">www.bioworlde.com</a>

Bioworld technology, co. Ltd. (China)



#### X. TYPICAL DATA

The standard curve is for demonstration only. A standard curve must be run with each assay.

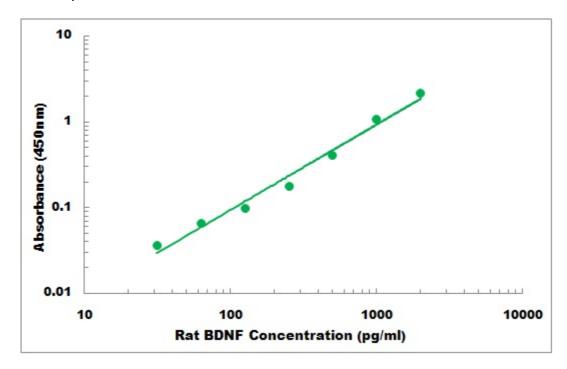

#### XI. SENSITIVITY

The minimum detectable dose of Rat BDNF is typically less than 2 pg/ml.

#### XII. SPECIFICITY

The Rat BDNF ELISA Kit allows for the detection and quantification of endogenous levels of natural and/or recombinant Rat BDNF proteins within the range of 31.2 pg/ml - 2000 pg/ml.

12

Bioworld Technology, Inc. (USA)

Email: <a href="mailto:info@bioworlde.com">info@bioworlde.com</a>
Web: <a href="mailto:www.bioworlde.com">www.bioworlde.com</a>

Bioworld technology, co. Ltd. (China)

XIII. **CROSS REACTIVITY** 

No detectable cross-reactivity with other relevant proteins.

XIV. REFERENCES

1. Jones, K. R.; Reichardt, L. F. Molecular cloning of a human gene that is a member

of the nerve growth factor family. Proc. Nat. Acad. Sci. 87: 8060-8064, 1990.

2. Kovalchuk, Y.; Hanse, E.; Kafitz, K. W.; Konnerth, A. Postsynaptic induction of

BDNF-mediated long-term potentiation. Science 295: 1729-1734, 2002.

3. Lyons, W. E.; Mamounas, L. A.; Ricaurte, G. A; Coppola, V.; Reid, S. W.; Bora, S. H.;

Wihler, C.; Koliatsos, V. E.; Tessarollo, L. Brain-derived neurotrophic factor-deficient

mice develop aggressiveness and hyperphagia in conjunction with brain serotonergic

abnormalities. Proc. Nat. Acad. Sci. 96: 15239-15244, 1999.

4. Rios, M.; Fan, G.; Fekete, C.; Kelly, J.; Bates, B.; Kuehn, R.; Lechan, R. M.; Jaenisch, R.

Conditional deletion of brain-derived neurotrophic factor in the postnatal brain leads

to obesity and hyperactivity. Molec. Endocr. 15: 1748-1757, 2001.

5. Baker-Herman, T. L.; Fuller, D. D.; Bavis, R. W.; Zabka, A. G.; Golder, F. J.; Doperalski,

N. J.; Johnson, R. A.; Watters, J. J.; Mitchell, G. S. BDNF is necessary and sufficient for

spinal respiratory plasticity following intermittent hypoxia. Nature Neurosci. 7: 48-55,

2004.

13

Bioworld Technology, Inc. (USA) Bioworld technology, co. Ltd. (China)

Email: info@bioworlde.com Email: info@biogot.com Web: www.bioworlde.com

Web: www.biogot.com

## XV. TROUBLESHOOTING GUIDE

Problem	Possible Cause	Solution
High signal and background in	Insufficient washing	Increase number of washes
all wells		Increase time of soaking
		between in wash
	Too much Streptavidin-HRP	Check dilution, titration
	Incubation time too long	Reduce incubation time
	Development time too long	Decrease the incubation
		time before the stop solution
		is added
No signal	Reagent added in incorrect	Review protocol
	order, or incorrectly prepared	
	Standard has gone bad (If	Check the condition of
	there is a signal in the sample	stored standard
	wells)	
	Assay was conducted from an	Reagents allows to come to
	incorrect starting point	20 - 30 °C before performing
		assay
Too much signal-whole plate	Insufficient washing-unbound	• Increase number of washes
turned uniformly blue	Streptavidin-HRP remaining	Carefully
	Too much Streptavidin-HRP	Check dilution
	Plate sealer or reservoir	Use fresh plate sealer and
	reused, resulting in presence of	reagent reservoir for each
	residual Streptavidin-HRP	step
Standard curve achieved but	Plate not developed long	Increase substrate solution
poor discrimination between	enough	incubation time
point	Improper calculation of	Check dilution, make new
	standard curve dilution	standard curve
No signal when a signal is	Sample matrix is masking	More diluted sample
expected, but standard curve	detection	Recommended
looks fine		
Samples are reading too high,	Samples contain protein levels	Dilute samples and run
but standard curve is fine	above assay range	Again
Edge effect	Uneven temperature around	Avoid incubating plate in
	work surface	areas where environmental
		conditions vary