

# PRODUCT DATA SHEET

Bioworld Technology,Inc.

# Cytokeratin 18 monoclonal antibody

Catalog: MB65912 Host: Mouse Reactivity: Human, Mouse, Rat

#### **BackGround:**

Keratins (cytokeratins) are intermediate filament proteins that are mainly expressed in epithelial cells. Keratin heterodimers composed of an acidic keratin (or type I keratin, keratins 9 to 23) and a basic keratin (or type II keratin, keratins 1 to 8) assemble to form filaments. Keratin isoforms demonstrate tissue- and differentiation-specific profiles that make them useful as research biomarkers. Research studies have shown that mutations in keratin genes are associated with skin disorders, liver and pancreatic diseases, and inflammatory intestinal diseases.

#### **Product:**

Liquid in 0.42% Potassium phosphate, 0.87% Sodium chloride, pH 7.3, 30% glycerol, and 0.01% sodium azide.

#### **Molecular Weight:**

~ 48 kDa

## **Swiss-Prot:**

P05783

### **Purification&Purity:**

The antibody was affinity-purified from mouse antiserum by affinity-chromatography using epitope-specific immunogen and the purity is > 95% (by SDS-PAGE).

## **Applications:**

WB (1/1000 - 1/3000), IF/ICC (1/100 - 1/200)


#### Storage&Stability:

Store at  $4\,\mathrm{C}$  short term. Aliquot and store at -20  $\mathrm{C}$  long term. Avoid freeze-thaw cycles.


## **Specificity:**

Recognizes endogenous levels of Cytokeratin 18 protein.

## DATA:

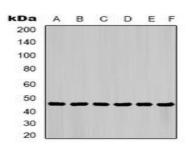

Western blot analysis of Cytokeratin 18 expression in HepG2 (A), Hela (B), mouse liver (C), mouse skeletal muscle (D), C2C12 (E), rat heart (F) whole cell lysates.

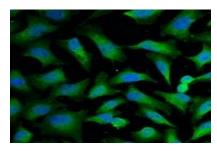

Immunofluorescent analysis of Cytokeratin 18 staining in Hela cells. Formalin-fixed cells were permeabilized with 0.1% Triton X-100 in TBS for 5-10 minutes and blocked with 3% BSA-PBS for 30 minutes at room temperature. Cells were probed with the primary antibody in 3% BSA-PBS and incubated overnight at 4 °C in a hidified chamber. Cells were washed with PBST and incubated with a FITC-conjugated secondary antibody (green) in PBS at room temperature in the dark. DAPI was used to stain the cell nuclei (blue).

## Note:

For research use only, not for use in diagnostic procedure.

# Bioworld Technology, Inc.

Add: 1660 South Highway 100, Suite 500 St. Louis Park,

MN 55416,USA.

Email: <u>info@bioworlde.com</u>

Tel: 6123263284 Fax: 6122933841 Bioworld technology, co. Ltd.

Add: No 9, weidi road Qixia District Nanjing, 210046,

P. R. China.

Email: <u>info@biogot.com</u>
Tel: 0086-025-68037686
Fax: 0086-025-68035151