

PRODUCT DATA SHEET

Bioworld Technology,Inc.

PIP monoclonal antibody

Catalog: MB66213 Host: Mouse Reactivity: Human

BackGround:

Gross cystic disease of the breast is a benign premenopausal disorder in which cysts are a predominant pathological lesion. These cysts appear to be formed from excessive aporcine cystic secretions. This fluid is composed of several glycoproteins including a unique 15kD monomer protein, PIP. Cytosolic analysis of normal tissue specimens from all major organs has demonstrated PIP in aporcine epithelia, lacrimal, ceruminous and Moll's glands and in numerous serous cells of the submandibular, tracheal, bronchial, sublingual and minor salivary glands. PIP and prostate specific antigen are co-expressed in androgen receptor-positive breast tumors.

Product:

Mouse IgG1. Liquid in PBS containing 50% glycerol, 0.2% BSA and 0.01% sodium azide.

Molecular Weight:

~ 16 kDa

Swiss-Prot:

P12273

Purification&Purity:

The antibody was affinity-purified from mouse antiserum by affinity-chromatography using epitope-specific immunogen and the purity is > 95% (by SDS-PAGE).

Applications:

IHC (1/100 - 1/300)


Storage&Stability:

Store at $4\,\mathrm{C}$ short term. Aliquot and store at $-20\,\mathrm{C}$ long term. Avoid freeze-thaw cycles.


Specificity:

Recognizes endogenous levels of PIP protein.


DATA:

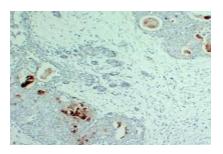

Immunohistochemical analysis of PIP staining in human breast carcinoma formalin fixed paraffin embedded tissue section. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH 6.0). The section was then incubated with the antibody at room temperature and detected using an HRP conjugated compact polymer system. DAB was used as the chromogen. The section was then counterstained with haematoxylin and mounted with DPX.

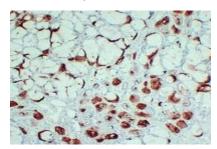

Immunohistochemical analysis of PIP staining in human salivary gland formalin fixed paraffin embedded tissue section. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH 6.0). The section was then incubated with the antibody at room temperature and detected using an HRP conjugated compact polymer system. DAB was used as the chromogen. The section was then counterstained with haematoxylin and mounted with DPX.

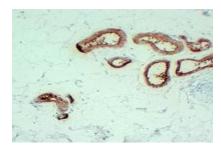

Immunohistochemical analysis of PIP staining in human skin formalin

Bioworld Technology, Inc.

Add: 1660 South Highway 100, Suite 500 St. Louis Park,

MN 55416,USA.

Email: <u>info@bioworlde.com</u>

Tel: 6123263284 Fax: 6122933841

Bioworld technology, co. Ltd.

Add: No 9, weidi road Qixia District Nanjing, 210046,

P. R. China.

Email: <u>info@biogot.com</u>
Tel: 0086-025-68037686
Fax: 0086-025-68035151



PRODUCT DATA SHEET

Bioworld Technology,Inc.

fixed paraffin embedded tissue section. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH 6.0). The section was then incubated with the antibody at room temperature and detected using an HRP conjugated compact polymer system. DAB

was used as the chromogen. The section was then counterstained with haematoxylin and mounted with DPX.

Note:

For research use only, not for use in diagnostic procedure.

Bioworld Technology, Inc.

Add: 1660 South Highway 100, Suite 500 St. Louis Park,

MN 55416,USA.

Email: info@bioworlde.com

Tel: 6123263284 Fax: 6122933841 Bioworld technology, co. Ltd.

Add: No 9, weidi road Qixia District Nanjing, 210046,

P. R. China.

Email: <u>info@biogot.com</u>
Tel: 0086-025-68037686
Fax: 0086-025-68035151