

# **Beta-tubulin Rabbit monoclonal antibody**

Catalog: MB66321

Host: Rabbit

Reactivity: Human, Mouse, Rat

## **BackGround:**

The cytoskeleton consists of three types of cytosolic fibers: microtubules, microfilaments (actin filaments), and intermediate filaments. Globular tubulin subunits comprise the microtubule building block, with  $\alpha/\beta$ -tubulin heterodimers forming the tubulin subunit common to all eukaryotic cells.  $\gamma$ -tubulin is required to nucleate polymerization of tubulin subunits to form microtubule polymers. Many cell movements are mediated by microtubule action, including the beating of cilia and flagella, cytoplasmic transport of membrane vesicles, chromosome alignment during meiosis/mitosis, and nerve-cell axon migration. These movements result from competitive microtubule polymerization and depolymerization or through the actions of microtubule motor proteins.

**Product:** 

Liquid in 50mM Tris-Glycine (pH 7.4), 0.15M NaCl, 50% Glycerol, 0.01% Sodium azide and 0.05% BSA.

**Molecular Weight:** 

~ 55 kDa

**Swiss-Prot:** 

P07437

**Purification&Purity:** 

The antibody was purified by immunogen affinity chromatography.

## **Applications:**

WB (1/500 - 1/1000), IHC (1/50 - 1/100), IP (1/10 - 1/50)


# Storage&Stability:

Store at  $4 \,^{\circ}{\rm C}$  short term. Aliquot and store at  $-20 \,^{\circ}{\rm C}$  long term. Avoid freeze-thaw cycles.


## **Specificity:**

Recognizes endogenous levels of Beta-tubulin protein.

#### **DATA:**

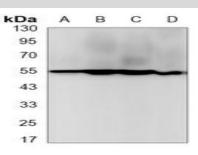

Western blot analysis of Beta-tubulin expression in rat brain (A), C6 (B), NIH3T3 (C), Hela (D) whole cell lysates.

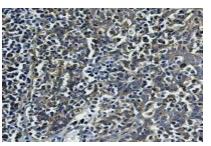

Immunohistochemical analysis of Beta-tubulin staining in human tonsil formalin fixed paraffin embedded tissue section. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH 6.119). The section was then incubated with the antibody at room temperature and detected using an HRP conjugated compact polymer system. DAB was used as the chromogen. The section was then counterstained with haematoxylin and mounted with DPX.

#### Note:

For research use only, not for use in diagnostic procedure.

#### Bioworld Technology, Inc.

 Add:
 1660 South Highway 100, Suite 500 St. Louis Park, MN 55416,USA.

 Email:
 info@bioworlde.com

 Tel:
 6123263284

 Fax:
 6122933841

# Bioworld technology, co. Ltd.

 
 Add:
 No 9, weidi road Qixia District Nanjing, 210046, P. R. China.

 Email:
 info@biogot.com

 Tel:
 0086-025-68037686

 Fax:
 0086-025-68035151