

## PRODUCT DATA SHEET

Bioworld Technology,Inc.

# **ALPP monoclonal antibody**

Catalog: MB66638 Host: Mouse Reactivity: Human

#### **BackGround:**

Mammalian alkaline phosphatases (APs) are highly conserved zinc-containing allosteric enzymes that are able to hydrolyze and transphosphorylate a wide range of compounds. There are four known human alkaline phosphaisozymes: **TNAP** (tissue-nonspecific; tase bone/liver/kidney), ALPP (placental), ALPP2 (germ cell), and ALPI (intestinal). Placental alkaline phosphatase (ALPP) is bound to the plasma membrane via a glycosyl-phosphatidylinositol (GPI) anchor. It is expressed primarily in the placenta and may be involved in transplacental IgG transport. ALPP has been found to be overexpressed on the surface of several different types of solid tumor cells and elevated serum concentrations of ALPP and ALPP-like enzymes has been found to be associated with ovarian, cervical, and testicular cancer.

## **Product:**

Liquid in PBS containing 50% glycerol, 0.5% BSA and 0.02% sodium azide, pH 7.3.

## **Molecular Weight:**

~ 70 kDa

### **Swiss-Prot:**

P05187

## **Purification&Purity:**

The antibody was purified by immunogen affinity chromatography.

## **Applications:**

WB (1/500 - 1/1000), IHC (1/50 - 1/100)


#### **Storage&Stability:**

Store at  $4\,\mathrm{C}$  short term. Aliquot and store at -20  $\mathrm{C}$  long term. Avoid freeze-thaw cycles.


#### **Specificity:**

Recognizes endogenous levels of ALPP protein.

#### **DATA:**

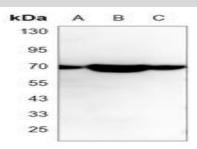

Western blot analysis of ALPP expression in HepG2 (A), A431 (B), MCF7 (C) whole cell lysates.

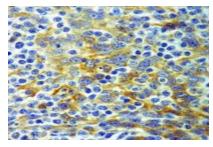

Immunohistochemical analysis of ALPP staining in human tonsil formalin fixed paraffin embedded tissue section. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH 6.141). The section was then incubated with the antibody at room temperature and detected using an HRP conjugated compact polymer system. DAB was used as the chromogen. The section was then counterstained with haematoxylin and mounted with DPX.

#### Note:

For research use only, not for use in diagnostic procedure.

## Bioworld Technology, Inc.

Add: 1660 South Highway 100, Suite 500 St. Louis Park,

MN 55416,USA.

Email: <u>info@bioworlde.com</u>

Tel: 6123263284 Fax: 6122933841 Bioworld technology, co. Ltd.

Add: No 9, weidi road Qixia District Nanjing, 210046,

P. R. China.

Email: <u>info@biogot.com</u>
Tel: 0086-025-68037686
Fax: 0086-025-68035151